

**labodam**



**PHARMACY REFRIGERATOR  
LPRF-A14**

# PHARMACY REFRIGERATOR LPRF-A14


Pharmacy Refrigerator LPRF-A14 is a digital microprocessor controlled temperature refrigerator with holding capacity of 656 liters and adjustable design storage shelves. Features Inbuilt full automatic defrost system with CFC-Free refrigerant. It has standard electrical heater with double glass observation door and Noiseless, user friendly operation mode. Auto audio-visual signal alarm protection and it temperature ranges from 2 to 8°C with the LED display.

## FEATURES

- ◆ Accurate temperature control to make inner temperature maintain 2 to 8 °C and increment at 0.1°C
- ◆ High-tech integrated sensors to display and control temperature
- ◆ Auto-defrost to remove moisture on cooling surface
- ◆ Removable outdoor double seal design for better temperature maintain and energy saving
- ◆ Microprocessor controlled forced-air cooling system
- ◆ Digital temperature displays with 0.1°C resolution
- ◆ Equipped with double layer door
- ◆ CFC-Free refrigerant and Shelves adjustable design
- ◆ Epoxy coated exterior anti-bacteria and easy cleaning
- ◆ Caster with break design for easy moving and lock
- ◆ Interior LED light for easy observation
- ◆ Electrical heater to decrease condensation on door in high humidity environment
- ◆ Equipped by Five alarm conditions High/Low Temperature; Power Failure; Sensor Error; Door Ajar, Low Voltage in Backup Battery.

## APPLICATIONS

The pharmacy refrigerator is suitable for installation in drug stores, pharmacy companies, hospitals, epidemic prevention centre's to store medicines, testing kits, reagents, vaccines etc.

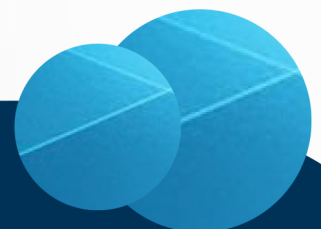

## SPECIFICATION

<b>Model</b>	<b>LPRF-A14</b>
Cabinet Type	Upright
Gross Capacity	656 L
Climate Class	10 to 32°C
Cooling Type	Forced Air Cooling
Defrost Mode	Automatic
Refrigerant	CFC-Free
Door	Double Glass Door
Noise	42 dB
Temperature Range	2 to 8°C
Uniformity	± 3°C
Temperature Control	Microprocessor
Display	LED
Internal Dimension	1100 × 454 × 1325 mm
External Dimension	1220 × 642 × 1885 mm
Packing Dimension	1290 × 695 × 1940 mm
Power Voltage	220 to 240 / 50 to 60 V/Hz
Power Supply	443 W
Net/Gross Weight	158 /188 kg

## STANDARD ACCESSORIES

<b>Accessories No.</b>	<b>Name</b>	<b>Qty</b>
1	Caster	4
2	Foot	Y
3	Shelf/Basket	10/-

